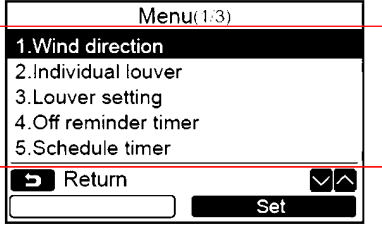
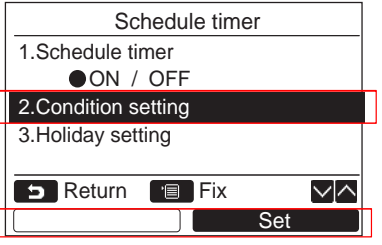
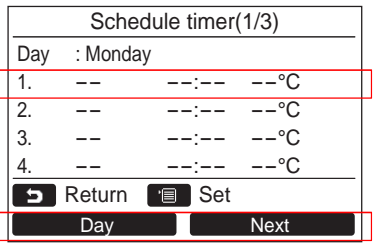
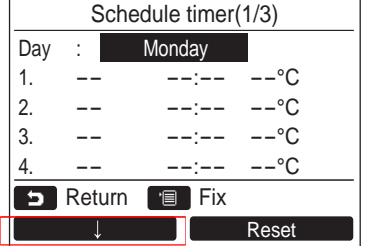
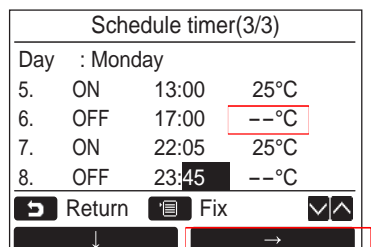
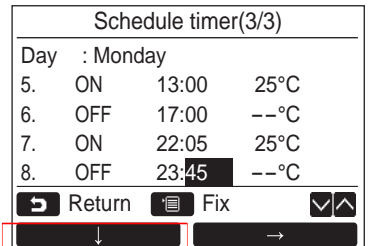


<p>Press the [MENU] button to display the "Menu" screen</p> <p>Press the [] [] buttons to scroll through settings & select option "5 Scheduled timer" then press F2 Set</p>	
<p>Press the [] [] buttons to scroll through settings to select option "2 Condition setting" then press F2 Set</p>	
<p>The current settings are displayed</p> <p>Press the F1 Day to confirm the day settings</p> <p>Press the F2 Next to confirm setting, 8 settings appear</p>	
<p>Press the [MENU] button</p> <p>Press the [] [] buttons to select the day to set</p> <p>Press the F1 button</p>	
<p>Press the [] [] buttons to select "ON" or "OFF"</p> <p>Select "ON" to set start time and set temperature settings</p> <p>Select "OFF" to set stop time.</p> <p>"- -" indicates that item has not been set</p> <p>Press the F2 button to select time or temperature</p> <p>If "- -" is displayed (ON/OFF not set) time or temperature cannot be set</p>	
<p>Press the [] [] buttons to set time or temperature</p> <p>Press the F1 to program next sequence. Up to 8 sequence settings per day can be programmed</p> <p>Press the [MENU] button to display day selection screen</p> <p>Press the [] [] buttons to select the next day to set</p> <p>Repeat procedures above to program day, time & temperature settings</p>	
<p>Press the [MENU] button</p> <p>Press the F1 Yes to confirm</p> <p>Press the F2 No to return to setting screen</p>	